
Extremely easy to use clinical gait analysis for rehab purposes









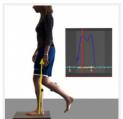

Highlights

Accurate data acquisition

Data is accurately acquired from high quality sensors

High Picture Quality

Through the use of High Quality camera's and video acquisition software high picture quality enables easy data analysis

Record and analysis in one software package

One comprehensive software package for the complete analysis: recording, analysing and reports

Customized reporting included

You can create your customized reports with your own logo and the data you want within the software.

Extendable configuration

The standard CMAX setup can be easily extended with extra camera's, EMG and forceplates





Extremely easy to use clinical gait analysis for rehab purposes



CMAX is a turnkey solution for 2 dimensional gait analyses. This complete solution for video gait analysis is especially developed to use in a clinical setting and contains hardware as well as software. The combination of the hardware (video, force plate and potential EMG) ensures that the required data for gait analysis can be collected, while the software combines all the data and contains some tools for easy analysis. Gait analysis is mainly used in rehabilitation and orthopedics to evaluate gait, determine the best treatment and evaluate the treatment.

The system is optimally designed for 2D video gait analysis of human lower extremities and trunk for e.g.

- children with cerebral palsy
- people with acquired brain impairment like CVA / stroke
- neuromuscular diseases like Polyomelites, Multiple sclerosis, Duchenne disease etc.
- fitting of ortheses and protheses

The system contains:

- workstation with large LCD screen, keyboard and mouse
- pre-installed software
- licenses
- cables
- frontal: first camera with lens and tripod or wallmount
- sagital: second camera with lens and tripod or wallmount
- forceplate
- calibration frame

The system can be easily extended and synchronised with an additional camera, EMG and an additional forceplate.

Features



Easy installation

Lode Biometry products are easy to install



Flexible Protocols

With Lode Biometry software products you can generate and execute customized protocols in a flexible way.



Automatic Synchronisation

All data streams are automatically synchronized for optimal data accuracy



Easy Calibration

Lode clinical gait analysis products can be easily calibrated with the calibration frame



Single Power Switch On

Through the special power unit all data acquisition devices can be switched with one press on a button



Moxie viewer compatible

The data is compatible with the free Moxie viewer 1 and 2, so it can be viewed on other locations.



Optimal room use

By smart lens selection smaller rooms can also be used for clinical gait analysis with CMAX



Extremely easy to use clinical gait analysis for rehab purposes



CMAX including Forceplate can a.o be extended with the following options:

Medical Keyboard

Maximum hygiene



Partnumber: 945860

Medical Mouse

Maximum hygiene



Partnumber: 945861

Additional Camera

Front and rear view in one recording



Partnumber: 990902

Additional Forceplate

Easy M-wave capturing



Partnumber: 990906

8 channel wireless EMG

Accurate measurement of timing and force of



Partnumber: 990903

16 channel wireless EMG

Accurate measurement of timing and force of



Partnumber: 990904

Walkway for CMAX

Easy installation without construction



Partnumber: 990908

CMAX rolling camera rail system

Optimal camera placement



Partnumber: 990825

Wall mount for Camera

Optimal camera placement



Partnumber: 990722



Extremely easy to use clinical gait analysis for rehab purposes



Specifications

| Accuracy | | |
|-------------------------------|---------|-----------|
| Sampling frequency | 1000 Hz | |
| Sensor surface length | 50 cm | 19.7 inch |
| Sensor surface width | 50 cm | 19.7 inch |
| Comfort | | |
| Remote Battery Indication | ✓ | |
| User Interface | | |
| English user interface | ✓ | |
| Customized Reports | ✓ | |
| Software available in English | ✓ | |
| Connectivity | | |
| USB connector | ✓ | |
| Dimensions | | |
| Walking surface length | 800 cm | 315 inch |
| Walking surface width | 120 cm | 47.2 inch |
| Product weight | 45 kg | 99.2 lbs |
| Compatibility | | |
| Microsoft Windows 10 | ✓ | |
| Microsoft Windows 8 | ✓ | |
| Microsoft Windows 7 | ✓ | |
| Included parts | | |
| PC included | ✓ | |
| PC software included | ✓ | |
| Windows 7 installed | ✓ | |
| Keyboard and mouse included | ✓ | |
| LCD screen included | ✓ | |
| All cables included | ✓ | |
| | | |

Order info

Partnumber: 990900



Lode B.V.

^{*}Specifications are subject to change without notice.